

CATALOGUE



YOUR COMPLETE VETERINARY SUPPLY AND SUPPORT SOLUTION



#### **ABOUT EUROVETS**

#### EUROVETS IS YOUR COMPLETE VETERINARY SUPPLY & SUPPORT SOLUTION

Eurovets, established in 2001 has become the Middle East's foremost supplier of some of the world's best known veterinary products to veterinary clinics throughout the UAE.

We represent leading multi-national brands in nutrition, veterinary equipment, medicines and consumables for companion, -equine and production animals. In addition to the traditional wholesaling, Eurovets provides warehousing, distribution, marketing and sales resources in order to support the product in the region.

We also offer a host of value-added services such as maintenance and support for all equipment supplied; service-based offerings such as Spectrum allergy testing; even end-to-end IT solutions for clinics.

Eurovets supports Continued Professional Development, imparting knowledge of the latest developments in the field to veterinary professionals through high-powered seminars.

We pride ourselves on our professional approach and offer vet specialized product advice where required.

In conclusion our mission is to not only deliver excellent products but to be committed to deliver exceptional customer service to all veterinary clinics and pet retailers in the region.

Eurovets YOUR headquarters for veterinary supply and support and a proud distributor of the IDEXX brand.




#### ABOUT IDEXX LABORATORIES

IDEXX Laboratories, Inc. develops, manufactures, and distributes various products and services for the companion animal veterinary and livestock markets worldwide. The company offers point-of-care veterinary diagnostic products comprising instruments and consumables, such as blood and urine chemistry, hematology, quantitative immunoassay testing, and urinalysis analyzers; and rapid assays.

The company was founded in 1983 and is headquartered in Westbrook, Maine and has been operating in the Middle east for over a decade.

IDEXX is the single most widely used Laboratory equipment amongst veterinarians in the region.

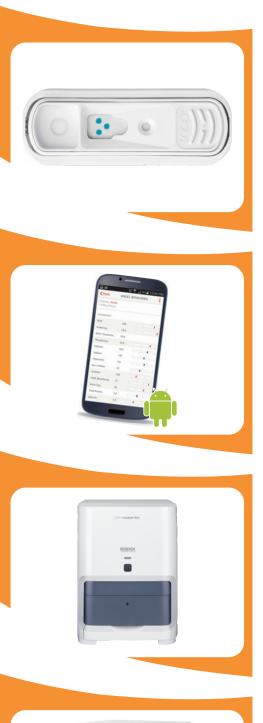

## INDEX

Page
Our favourites
The IDEXX VetConnect® PLUS mobile app
The SNAP Pro <sup>®</sup> Mobile Device
The Catalyst One <sup>™</sup> Chemistry Analyzer
IDEXX In-house Analyzers 6
Diagnostic Tests for IDEXX In-house Analyzers 8
IDEXX SNAP® Tests
IDEXX Education and Training
STEPs – Simple Tools for Effective Protocols





# OUR FAVOURITES

## Personalized alerts for a comprehensive and timely diagnostic picture

#### The IDEXX VetConnect® PLUS mobile app

- Access your diagnostics from IDEXX in-house analyzers, SNAP<sup>®</sup> tests and IDEXX Reference Laboratories anywhere, anytime.
- Receive a notification the moment diagnostic results are available by leaving Push Notifications on.
- Apply personalization so each veterinarian will only receive results for their own patients.

#### Download the app today.

Android device? Visit idexx.com/vcpmobile1 or the Google Play store. iPhone device? Visit idexx.com/vcpmobile or the Apple App Store.



## Increase the value that your SNAP<sup>®</sup> tests offer patients, clients and your practice

### The SNAP Pro® Mobile Device

Streamline your technician's work flow with the device that automatically activates SNAP tests, displays clear test results and sends them to the patient record – all at the touch of a button.

- Simple design and easy-to-read display improves client engagement, communication and compliance.
- Electronically captures SNAP test invoice charges to reduce lost revenue and improve profits.
- Instantly update the patient record for customized, client-friendly reports with IDEXX VetConnect® PLUS.



## A comprehensive workup at the point of care with real-time results

#### The Catalyst One<sup>™</sup> Chemistry Analyzer

- In one run, with one whole blood sample, get chemistry and electrolyte results in 8 minutes plus total  $T_4$  results in as little as 15.
- Expanding 30-test menu includes advanced diagnostics, like fructosamine, phenobarbital and Urine P:C Ratio for more insight.



### **IDEXX In-house Analyzers**

Deliver real-time care and increase practice revenue with a flexible, broad range of in-house profiles, a more efficient work flow and better information management.

### **IDEXX SmartService<sup>™</sup> Solutions**

IDEXX SmartService<sup>™</sup> Solutions is the highest level of customer support and service in the veterinary industry. With electronic troubleshooting, proactive maintenance and automatic software upgrades, it helps you increase your instrument uptime and improve practice efficiency.

#### **IDEXX** integration technology

IDEXX integration technology allows data to flow back and forth between your practice management software and IDEXX in-house analyzers. This "two-way conversation" streamlines work flow by consolidating patient diagnostic data, minimizing data entry and typing errors and capturing lost charges. Activate IDEXX VetConnect<sup>®</sup> PLUS and view all your diagnostic results anytime, anywhere.

IDEXX VetLab <sup>®</sup> Station Laboratory Information Management System (LIMS)	SNAP Pro <sup>®</sup> Mobile Device	Catalyst Dx◎ Chemistry Analyzer	Catalyst One™ Chemistry Analyzer	VetStat◎ Electrolyte and Blood Gas Analyzer
<ul> <li>INFORMATION</li> <li>Increase your laboratory information management and diagnostic capabilities.</li> <li>User-friendly touch screen simplifies laboratory work flow</li> <li>Integration with most practice management software</li> <li>All patient results in a single easy-to-read report</li> <li>Comprehensive patient records and diagnostic history, including parameter trending</li> <li>Access to analyte descriptions and IDEXX VetLab Station snippet videos through the online Help system</li> <li>Remote instrument troubleshooting with IDEXX SmartService™ Solutions</li> </ul>	<ul> <li>MANAGEMENT</li> <li>Manage your SNAP* results and improve work flow and profits when running the SNAP tests you rely on every day.</li> <li>Automatically activates test device and displays results for interpretation to streamline technician's work flow</li> <li>Sends SNAP test results to the patient record at the touch of a button</li> <li>Electronically captures SNAP test invoice charges to reduce lost revenue when used with select practice management software</li> <li>Shares information with your practice management system</li> </ul>	<ul> <li>Revolutionary speed, throughput and flexibility put results in your hands more quickly so you can manage more cases simply and efficiently. Trusted IDEXX dry-slide technology offers accurate results.</li> <li>Preloaded profiles in easy-to-use CLIPs</li> <li>Onboard whole blood separation</li> <li>Multiple-patient load-and-go results</li> <li>Complete menu flexibility: 35 parameters (29 individual tests and 6 calculations) and 7 convenient preloaded CLIPs</li> </ul>	<ul> <li>CHEMISTRY</li> <li>Deliver personalized, exceptional care with core diagnostic results in one sample run.</li> <li>Share comprehensive results (total T4, chemistry and electrolytes) during the patient visit</li> <li>Get real-time answers for preanesthetic screens as well as sick- and well-patient visits</li> <li>Run just the tests you need with complete single-slide testing flexibility</li> <li>Save technician time with load-and-go efficiency</li> <li>See all of a patient's diagnostic results on one screen with IDEXX VetConnect<sup>®</sup> PLUS</li> </ul>	<ul> <li>Assess electrolytes, blood gases, acid-base balance, ionized calcium, glucose, anion gap, tCO<sub>2</sub> and bicarbonate in-house.</li> <li>Designed and validated specifically for veterinary use</li> <li>Convenient cassette choices to run plasma, serum, venous whole blood or arterial whole blood</li> <li>Rechargeable battery for up to four hours of independent operation</li> </ul>

## See all your diagnostic results on a single screen, anytime.

Only IDEXX gives you one comprehensive view of patient diagnostic results from all your in-house and reference laboratory testing:

- Consolidated historical results and trending graphs with IDEXX VetConnect® PLUS
- Comprehensive point-of-care chemistry and hematology solutions with real-time results
- Broadest menu of advanced reference laboratory diagnostics with fast turnaround times
- Automatic alerts with personalized notifications are sent to iPhone<sup>®</sup> or Android<sup>™</sup> devices the instant results are ready



				e l'
SNAPshot Dx <sup>®</sup> Analyzer	ProCyte Dx <sup>®</sup> Hematology Analyzer	LaserCyte◎ Dx Hematology Analyzer	Coag Dx™ Analyzer	IDEXX VetLab◎ UA™ Analyzer
<ul> <li>IMMUNOASSAY</li> <li>Add this analyzer to your in-house laboratory for immediate results and efficient management of a wide range of diseases. Trusted SNAP® ELISA technology.</li> <li>Reads multiple tests simultaneously and efficiently</li> <li>Use SNAP® Total T4, Cortisol and Bile Acids tests to help evaluate thyroid, adrenal and liver functions</li> <li>Now run SNAP® fPL™, SNAP® cPL™, SNAP® Feline Triple®, SNAP® FlV/FeLV Combo, SNAP® 4Dx® Plus, SNAP® Heartworm RT and SNAP® Feline proBNP tests</li> <li>Quickly screen and monitor dosages and effectiveness of treatments</li> </ul>	<ul> <li>Three cutting-edge technologies provide fast, accurate and complete results in just two minutes.</li> <li>An advanced five-part white blood cell differential with laser flow cytometry</li> <li>An absolute reticulocyte count is provided and feline platelets are analyzed using optical fluorescence</li> <li>Fast, precise count of red blood cells with Laminar Flow Impedance™</li> <li>Results validated for canine, feline, equine, bovine, ferret, rabbit, gerbil, guinea pig, pig and mini pig</li> <li>26 different parameters reported*, including band neutrophil and nRBC indicators</li> <li>Fluid analysis capabilities</li> </ul>	<ul> <li>HEMATOLOGY</li> <li>Laser flow cytometry providing in-depth, comprehensive CBC with a full five-part optical differential and an absolute reticulocyte count.</li> <li>Laboratory-quality laser flow cytometry technology individually examines each cell for accurate results</li> <li>Absolute reticulocyte counts are the gold standard for classifying anemia</li> <li>No daily maintenance</li> <li>Results validated for canine, feline, equine, ferret, pig and mini pig</li> <li>Fluid analysis capabilities</li> </ul>	<ul> <li>Detect and diagnose clotting disorders in seconds. Connects to the IDEXX VetLab Station for integrated results, electronic record keeping and parameter trending. Results are validated for veterinary use.</li> <li>Provides in-house testing for PT and aPTT in seconds for fresh or citrated whole blood</li> <li>LED optical detectors monitor blood flow to detect clot formation</li> <li>Lithium ion battery runs more than two hours on a full charge for portability</li> </ul>	<ul> <li>URINALYSIS</li> <li>Improves work flow and offers consistent results. Automatically reads IDEXX UA™ Strips and prints your results so there's no need to decipher strip colors.</li> <li>Simple operation: load strip, press a button, walk away—no need to wait for results</li> <li>Veterinary-specific strips validated for canine, feline and equine use</li> <li>Reflectance photometry technology uses LEDs to read test pads</li> </ul>

\*Number of reported parameters may vary according to species.

## **Diagnostic Tests for IDEXX In-house Analyzers**

Make confident real-time decisions. These fast, accurate and flexible veterinary-specific diagnostic tests give you comprehensive in-house capability.



Catalyst Dx<sup>®</sup> Chemistry Analyzer and Catalyst One<sup>™</sup> Chemistry Analyzer

		¢	A	¢	
Chem 17 CLIP	Chem 15 CLIP	Chem 10 CLIP	Lyte 4 CLIP	NSAID 6 CLIP	Equine 15 CLIP
The ideal comprehensive panel for canine testing (wellness, sick, emergency). Also for adult and geriatric canine preanesthetic screening. • 17 results per CLIP • 12 CLIPs per box ALB, ALB/GLOB, ALKP, ALT,		Screen young patients for safer anesthesia. • 10 results per CLIP • 12 CLIPs per box ALB, ALB/GLOB, ALKP, ALT, BUN, BUN/CREA, CREA, GLOB, GLU, TP	Easily add electrolytes to any analysis with the preloaded electrolyte CLIP. • 4 results per CLIP • 12 CLIPs per box CI, K, Na, Na/K	Establish a baseline, test during the trial period and continue monitoring your patients on NSAID therapy so you can adjust medication if necessary. • 6 results per CLIP • 12 CLIPs per box ALKP, ALT, AST, BUN,	Comprehensive equine testing in one convenient preloaded CLIP. • 15 results per CLIP • 12 CLIPs per box ALB, ALB/GLOB, ALKP, AST, BUN, BUN/CREA, Ca, CK, CREA, GGT, GLOB, GLU, LDH, TBIL, TP
AMYL, BUN, BUN/CREA, Ca, CHOL, CREA, GGT, GLOB, GLU, LIPA, PHOS, TBIL, TP	ALT, BUN, BUN/CREA, Ca, CHOL, CREA, GGT, GLOB,		Cannot be run as single slides: Cl, K, Na	BUN/CREA, CREA	

### 30 different tests with the Catalyst Dx<sup>®</sup> Chemistry Analyzer

Easily use single slides, create customized panels or use convenient preloaded CLIPs.

Alanine aminotransferase (ALT)
Albumin (ALB)
Alkaline phosphatase (ALKP)
Ammonia (NH3)
Amylase (AMYL)
Aspartate aminotransferase (AST)
Blood urea nitrogen (BUN)
Calcium (Ca)
Chloride (Cl)
Cholesterol (CHOL)

Creatine kinase (CK) Creatinine (CREA) Fructosamine (FRU) Gamma-glutamyltransferase (GGT) Glucose (GLU) Inorganic phosphate (PHOS) Lactate dehydrogenase (LDH) Lactate (LAC) Lipase (LIPA) Magnesium (Mg) Phenobarbital (PHBR) Potassium (K) Sodium (Na) Total bilirubin (TBIL) Total protein (TP) Total T4 (TT4) Triglycerides (TRIG) Uric acid (URIC) Urine creatinine (UCRE) Urine protein (UPRO)

Calculated values: Globulin (GLOB) Albumin/globulin (ALB/GLOB) ratio BUN/Creatinine (BUN/CREA) ratio Sodium/potassium (Na/K) ratio Urine protein:creatinine (UPC) ratio Osmolality


	•			Ŵ
Urine P:C Ratio	Fructosamine Slide	Phenobarbital Slide	Single Slides	QC CLIP
Detect and quantify proteinuria and diagnose early renal disease (dilution required; order diluent separately). • 2 chemistries per ratio • 6 ratios per box UCRE, UPRO, UPC	Confidently, quickly and cost-effectively monitor and manage diabetes in both cats and dogs. One sample provides rapid and accurate fructosamine levels. • 6 tests per box	<ul> <li>Monitor phenobarbital levels in your patients in minutes so you can treat with confidence.</li> <li>6 tests per box</li> <li>6 units of slide wash per box</li> </ul>	The ultimate in testing flexibility; run only the tests you want to run or create your own custom panels. • 12 tests per box ALB, ALKP, ALT, AMYL, AST, BUN, Ca, CHOL, CK, CREA, GGT, GLU, LAC, LDH, LIPA, Mg, NH3, PHOS, TBIL, TP, TRIC, URIC	of test results following the analyzer's monthly cleaning. • 6 CLIPs per box ALB, ALKP, ALT, Ca, GLU,
Cannot be run as single slides: UCRE, UPRO			trig, uric	

## Diagnostic Tests for IDEXX In-house Analyzers (continued)



VetStat<sup>®</sup> Electrolyte and Blood Gas Analyzer Cassettes

Electrolytes	Electrolytes 8 Plus	Respiratory/	Ionized Calcium	Glucose
Tests: Na⁺, K⁺, CI⁻	Tests: Na⁺, K⁺, CF, pH, <i>P</i> CO <sub>2</sub> , tCO <sub>2</sub> , HCO3 <sup>-</sup> and anion gap	Blood Gases Tests: Na <sup>+</sup> , K <sup>+</sup> , Cl <sup>-</sup> , pH, PC0 <sub>2</sub> , PO <sub>2</sub> , SO <sub>2</sub> , tHb,	Tests: Ca++	Tests: GLU
Use this cassette every time you run a biochemical profile to obtain a complete picture of your patient's health. • Get rapid, accurate results for major electrolytes • Uncover specific conditions, such as hyperkalemia and Addison's disease	<ul> <li>Obtain a complete metabolic/respiratory picture so you can individualize fluid therapy when needed and identify and treat the causes of acid-base disorders.</li> <li>Get rapid, accurate results for major electrolytes, acid-base and ventilation</li> <li>Uncover specific conditions, such as hyperkalemia and Addison's disease</li> <li>Monitor patients on therapeutic medications and fluid therapy</li> <li>Evaluate blood pH and determine if underlying causes are respiratory (airway obstruction, lung disease, central nervous system disorders) or metabolic (renal disease, lactic acidosis, diabetic ketoacidosis, dehydration, ethylene glycol toxicity)</li> </ul>	<ul> <li>PG2, PG2, SG2, HB, tCO2, HCO3<sup>-</sup>, base excess and anion gap</li> <li>Use this cassette for patients with respiratory and cardiovascular disorders, advanced anesthetic cases and all patients on ventilators to improve their chances of success in critical situations.</li> <li>Assess respiratory and breathing abnormalities</li> <li>Monitor cardiovascular and respiratory functions</li> <li>Assess a critical patient's need for supplemental oxygen immediately</li> <li>Monitor oxygen and respiratory functions</li> </ul>	<ul> <li>Confirm calcium</li> <li>imbalances on any patient</li> <li>with an abnormal total</li> <li>calcium, so you can begin</li> <li>treating your patients –</li> <li>and easing your clients'</li> <li>minds – as soon as</li> <li>possible.</li> <li>Measure the biologically</li> <li>active form of calcium</li> <li>in your patients' blood</li> <li>Assess increases or</li> <li>decreases in ionized</li> <li>calcium levels,</li> <li>which can indicate</li> <li>devastating disorders,</li> <li>such as lymphoma,</li> <li>acute pancreatitis,</li> <li>renal disease or an</li> <li>underactive parathyroid</li> <li>gland</li> <li>Monitor pregnant and</li> <li>lactating patients</li> </ul>	Quick assessments of glucose values will help you treat hypoglycemic and diabetic patients, as well as monitor patients on insulin therapy and oral hypoglycemics for diabetes. • Use whole blood to quickly assess glucose values

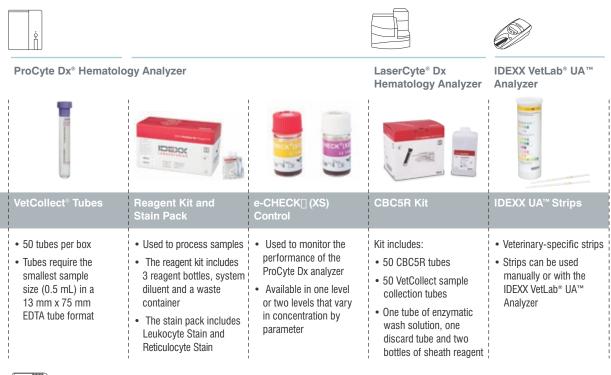

Coag Dx<sup>™</sup> Analyzer Cartridges

aPTT	Citrate aPTT	PT	Citrate PT
<ul> <li>Activated partial thromboplastin time test to screen for intrinsic and common coagulation pathway disorders.</li> <li>For use with fresh whole blood</li> <li>Smaller sample size; as little as 0.2 mL needed</li> <li>Must be run immediately</li> <li>Separate blood draw needed for each cartridge</li> </ul>	<ul> <li>Activated partial thromboplastin time test to screen for intrinsic and common coagulation pathway disorders.</li> <li>Perform test when convenient; citrated whole blood keeps at room temperature for up to two hours</li> <li>Test both citrate PT and citrate aPTT with one sample draw</li> <li>Use a 3.2% sodium citrate evacuated tube or syringe</li> </ul>	<ul> <li>Prothrombin time test to screen for extrinsic and common coagulation pathway disorders.</li> <li>For use with fresh whole blood</li> <li>Smaller sample size; as little as 0.2 mL needed</li> <li>Must be run immediately</li> <li>Separate blood draw needed for each cartridge</li> </ul>	<ul> <li>Prothrombin time test to screen for extrinsic and common coagulation pathway disorders.</li> <li>Perform test when convenient; citrated whole blood keeps at room temperature for up to two hours</li> <li>Test both citrate PT and citrate aPTT with one sample draw</li> <li>Use a 3.2% sodium citrate evacuated tube or syringe</li> </ul>

\*e-CHECK is a trademark or registered trademark of Sysmex Corporation.

## **IDEXX SNAP<sup>®</sup> Tests**


Point-of-care testing that provides the information you need to send more clients home with answers and the confidence that you're doing everything to care for their pets.

SNAP <sup>®</sup> Pro Mobile Devi	ce and SNAPshot Dx <sup>®</sup> A	nalyzer		SNAP <sup>®</sup> Pro Mobile Dev	ice
Case M	Carre M	Con. In	Con. In	Car. m	Con m
SNAP <sup>®</sup> 4Dx <sup>®</sup> Plus Test	SNAP <sup>®</sup> Heartworm RT Test	SNAP <sup>®</sup> Feline Triple <sup>®</sup> Test	SNAP <sup>®</sup> FIV/FeLV Combo Test	SNAP <sup>®</sup> Giardia Test	SNAP <sup>®</sup> Parvo Test
25°C 8 2°C 30	25°C 8 2°C 15 30	8°C 10 2°C 15 30		8 2°C 15	25°C 8 2°C 5
Provide a more thorough annual parasitic disease screen and reinforce the importance of tick control. The SNAP 4Dx Plus Test provides accurate results for 6 vector-borne disease—heartworm, Lyme disease, <i>Ehrlichia</i> <i>ewingii</i> , <i>Ehrlichia</i> <i>ewingii</i> , <i>Ehrlichia</i> <i>canis</i> , <i>Anaplasma</i> <i>phagocytophilum</i> and <i>Anaplasma platys</i> —in 8 minutes using one blood sample. <b>Sample requirement:</b> Serum, plasma or whole blood	Use this sensitive test for detecting even low worm burdens in just eight minutes. Minimal sample volume required. <b>Sample requirement:</b> Serum, plasma or whole blood	In communities where dogs get heartworm, cats can too. Raise the standard of care for cats in your practice with this three-in-one test. Features FIV and FeLV detection as well as a screen for feline heartworm antigen for the same price as our SNAP FIV/FeLV Combo Test. <b>Sample requirement:</b> Serum, plasma or whole blood	American Association of Feline Practitioners (AAFP) guidelines recommend testing all at-risk cats, sick cats and kittens to determine infection status. Feline immunodeficiency virus (FIV) and feline leukemia virus (FeLV) infections are two main causes of illness and death in cats. Screen for both diseases with just one blood sample. <b>Sample requirement:</b> Serum, plasma or whole blood	<i>Giardia</i> infection is often missed with fecal floats alone. Companion Animal Parasite Council (CAPC) recommends a combination of fecal flotation with centrifugation, direct smear and, for symptomatic patients, a <i>Giardia</i> ELISA. Catch more <i>Giardia</i> -positives with the quick and accurate SNAP <i>Giardia</i> Test. <b>Sample requirement:</b> Feces	Detect parvovirus during viral shedding before clinical symptoms appear. The SNAP Parvo Test will not cross-react with modified live vaccines, and easy-to-use swab greatly reduces risk of contamination. <b>Sample requirement:</b> Feces

Species. O Development time in minutes. Store at room temperature or in refrigerator. Store in refrigerator. Number of tests per kit. ALL components MUST be at room temperature before running the test.

KEY



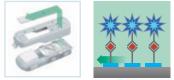

### **Trusted SNAP® ELISA technology**

The ELISA (enzyme-linked immunosorbent assay) is one of the most sensitive and reproducible diagnostic technologies available. IDEXX puts this superior technology to work in the small, easy-to-use format you've come to know as a SNAP test.

For detailed specifications on any SNAP test, visit idexx.com/snap. See a SNAP test animation at idexx.com/idexxelisa.

50



#### SNAP® Pro Mobile Device and SNAPshot Dx® Analyzer

ł

5

La.

SNAPshot Dx® Analyzer

#### ł (Contraction)

Con M	Con. M	Con. M	Con. M	Con. M	Con. M
SNAP <sup>®</sup> Feline proBNP Test	SNAP <sup>®</sup> cPL <sup>™</sup> Test (canine pancreas-specific lipase)	SNAP <sup>®</sup> fPL <sup>™</sup> Test (feline pancreas-specific lipase)	SNAP <sup>®</sup> Total T <sub>4</sub> Test[] SNAP <sup>®</sup> T <sub>4</sub> Test <sup>†</sup>	SNAP <sup>®</sup> Cortisol Test	SNAP <sup>®</sup> Bile Acids Test
	8°C 10 2°C 10	<b>1</b> 0 2 <b>1</b> 0 <b>21</b> 0	7 2°€ 6 15		7 2°C 6
The first pet-side cardiac biomarker blood test in the market, the SNAP Feline proBNP Test simply and quickly measures the level of NTproBNP that has been released by the heart. The level of NTproBNP in circulation provides insight into the amount of stress and stretch the heart is experiencing. <b>Sample requirement:</b> Serum or EDTA plasma	The only pet-side test to measure canine lipase levels specific to the pancreas, making it the most accurate point-of- care test to help diagnose or rule out pancreatitis in dogs presenting with vomiting, anorexia or abdominal pain. <b>Sample requirement:</b> Serum	The only pet-side test to measure feline lipase levels specific to the pancreas, making it the most accurate point-of- care test to help diagnose or rule out pancreatitis in cats presenting with lethargy, decreased appetite, dehydration, weight loss, vomiting and/or diarrhea. <b>Sample requirement:</b> Serum	Go from sample collection to results in just 15 minutes (not including incubation time). Screen for hyperthyroidism and hypothyroidism in one test, and begin treatment or follow-up testing in the same visit. Also, quickly monitor medications. <b>Sample requirement:</b> Serum or plasma The SNAP Total I Test is only compatible with the SNAPshot Dx Analyzer. 'The SNAP T4 Test is only compatible with the SNAP Reader.	Screen for Cushing's syndrome and Addison's disease and begin treatment in the same visit. Quantitative results make it easy to evaluate the effectiveness of treatments and adjust medications as needed. <b>Sample requirement:</b> Serum	The first fully quantitative in-house ELISA that helps you assess liver function and begin necessary treatment immediately. Quickly identify occult liver disease and detect possible portosystemic shunts, or use regularly to monitor response to treatment. <b>Sample requirement:</b> Serum

5

### **IDEXX Education and Training**

More than 50,000 veterinary practitioners, technicians and practice management professionals visit the IDEXX Learning Center every month for the most current veterinary information, advances in technology and medicine, and continuing education credits. Log in for knowledge you can put into practice today.



Log in to your account or register at idexxlearningcenter.com.

## **STEPs – Simple Tools for Effective Protocols**

Train your whole staff and build a stronger practice with STEPs. Whether you're learning to use a new instrument or implementing a new testing protocol, you'll find guides and informative resources to get everyone in your practice on the same track.

### Choose from these protocols:



Explore our growing list of protocols and other training resources at idexxlearningcenter.com/steps.

C # 3 trace data and an and a second place of the second sec	Home / Lab Draw / Tana		* 0 0 M 4	-
The State of the S		Contral Contral		
Max         Max         Max         Max           No         N         N         N         N	1 103 403 1 103 103 1 103 103 1 103 103	100 - 100 100 - 100 100 - 100		1
Alternational of Balling Constraining of Balling Const				-
Anno Anno Anno Anno Anno Anno Anno Anno				
Table In case of the second			1000	

### Take in the view of a greater vision of patient care



#### See all your diagnostic results on a single screen, anytime

Only IDEXX gives you one comprehensive view of patient diagnostic results from all your in-house and reference laboratory testing:

- Consolidated historical results and trending graphs with IDEXX VetConnect<sup>®</sup> PLUS
- Comprehensive point-of-care chemistry and hematology solutions with real-time results
- Broadest menu of advanced reference laboratory diagnostics with fast turnaround times
- Automatic alerts with personalized notifications are sent to iPhone<sup>®</sup> or Android<sup>™</sup> devices the instant results are ready

Information in these materials is subject to change without notice. No part of this document may be reproduced or transmitted in any form or by any means, electronic, mechanical or otherwise, for any purpose, without the express written permission of IDEXX Laboratories, Inc. IDEXX Laboratories may have patents or pending patent applications, trademarks, copyrights or other intellectual or industrial property rights covering this document or subject matter in this document. The furnishing of this document does not give a license to these property rights except as expressly provided in any written license agreement from IDEXX Laboratories.

PCR testing is a service performed pursuant to an agreement with Roche Molecular Systems, Inc.

Aller-g-complete is a registered trademark of GREER Laboratories, Inc. IPhone is a trademark of Apple Inc., registered in the United States and other countries. Android is a trademark of Google Inc. FedEx is a registered trademark of FedEx Corp. MWI is a registered trademark of MWI Veterinary Supply, Inc. All other %/TM marks are owned by IDEXX Laboratories, Inc. or its affiliates in the United States and/or other countries.

© 2014 IDEXX Laboratories, Inc. All rights reserved. • 09-67069-07 • The IDEXX Privacy Policy is available at idexx.com.

#### Talk with your IDEXX veterinary diagnostic consultant today to learn more.

#### ALL THE BEST BRANDS

Eurovets is dedicated to providing your clinic with all your veterinary needs. We offer a comprehensive product range at highly competitive prices for small and large animal practices, stables and wholesale outlets. This include biological, diagnostics, nutritionals, medicines, equipment, consumables and more.

#### Contact us today

#### **Eurovets Veterinary Suppliers**

Tel: +971 4 434 2436 **OR** Technical Support Team contact: +971 56 880 9212

E-mail: orders@eurovetsworld.com

Address: Warehouse A6, Dubai Science Park, Dubailand, Al Barsha South,

Dubai, United Arab Emirates



#### www.eurovetsworld.com

